

P-2000 Digital Polarimeter

Ultimate flexibility for a wide range of applications

*A new generation of polarimeters
offering simplicity, reliability and
complete validation protocols*



**tecnología
aplicada.com**

Jasco

P-2000 Digital Polarimeter

Ultimate flexibility for a wide range of applications ranging from quality control to research and technology



Foods

Pharmaceuticals

Sugars and Sweeteners

Essential oils, Flavors and Fragrances

Chemicals

JASCO, a leading manufacturer of polarimeters since 1967, is proud to introduce the new P-2000 multi-option polarimeter. The P-2000 is designed as a customizable polarimeter with various options for a range of applications and budgetary requirements. The instrument system can also be field upgraded as application requirements change. Options such as polarizers, wavelength filters, lamps and photomultiplier detectors provide a wide range of analytical wavelengths from UV-Vis to NIR.

Advanced instrument control

Two graphical user interfaces are available including a newly redesigned intelligent remote module (iRM) with a color LCD touch screen and Spectra Manager™ II software, the latest version of JASCO's innovative cross-platform spectroscopy software. Both of these control and analysis interfaces allow full-system control and advanced data processing. 21 CFR part 11 compliant versions for both interfaces can be ordered as an option.



The iRM conveniently guides the operator through routines from data acquisition to data processing. The obtained data can be automatically printed to USB printers, or saved to a compact flash memory card for further processing on a PC.

High speed, precision and accuracy

The P-2000 offers a response speed as fast as six degrees per second to provide reproducible data for each measurement with a resolution of 0.0001°. A wide dynamic range of up to ±90° enables the system to measure chiral compounds over a broad range of concentrations. The sample chamber is equipped with an integrated temperature probe; the current temperature is constantly updated on the instrument control screen.

Interference filter

A wide range of interference filters from UV-Vis to NIR

Light source switching mirror option

Light source

Up to two light sources can be installed. Available light sources are;

- WI (Tungsten-Halogen lamp)
- Na (Sodium lamp)
- Hg (Mercury lamp)


Polarizer

Faraday cell

I.Q. Accessory recognition

Analyzer

PMT detector

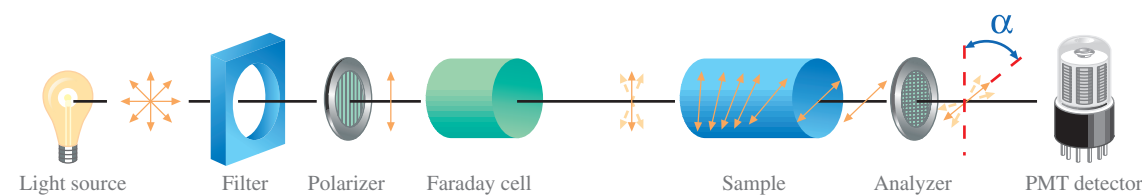

Full line of accessories

A full line of accessories, including a variety of cylindrical cells, a demountable stainless steel cell, Peltier thermostatted cell holder and monochromator options, are offered for functional expansion of the P-2000 to support a wide variety of applications. Quick-connect fittings with flow control valves in the sample compartment provide simple and secure exchange of water-thermostatted cells for temperature control with water baths. A newly developed Peltier cell holder offers accurate temperature control with ±0.1° accuracy.

Reliability and instrument verification

An automated validation program is available for GLP/cGMP laboratories or those regulated by the FDA. The program can be used for automated periodic validation of light source energy, zero repeatability, and rotational accuracy/repeatability. Inspected filters, sample cells with certified optical path length, and NIST-traceable rotation plates are also available.

Ultimate flexibility

With its single optical platform, the P-2000 can be customized to provide the optimum performance for every application. By selecting the most suitable combination of optical elements, the instrument can be used in the ultra-violet, visible and/or NIR regions.

Examples of versatile system configurations

A single wavelength system

546 nm with Hg or WI lamp

589 nm with Na or WI lamp

A dual wavelength system

589/546 nm with Na/Hg or WI lamp

5 wavelength system

589/578/546/436/365 nm with Na/Hg lamps

6 wavelength system

633/589/546/436/405/365 nm with WI lamp

7 wavelength system

633/589/578/546/436/405/365 nm with WI lamp

System expansion to UV range

334/313/302/296/280/254 nm with Hg lamp with the wavelength expansion option for the UV range

334/325 nm with WI lamp

with the wavelength expansion option for the UV range

Multiple wavelength system for UV/Vis range

633/589/578/546/436/405/365/334/325 nm with WI lamp with the wavelength expansion option for the UV range

System expansion to NIR range

880 nm with WI lamp with the wavelength expansion option for the NIR range

Multiple wavelength system for NIR and Visible range

880/589 nm with WI lamp with the wavelength expansion option for the NIR range

880/589/546 nm with WI lamp with the wavelength expansion option for the NIR range

Lamps and filters

Filter	Light Source		
	WI (Tungsten-Halogen lamp)	Na (Sodium lamp)	Hg (Mercury lamp)
880 mm	Quartz Faraday cell and R928		
633 mm			
589 mm			
578 mm			
546 mm			
436 mm			
405 mm			
365 mm			
334 mm	Quartz Faraday cell		Quartz Faraday cell
325 mm	Quartz Faraday cell		
313 mm			Quartz Faraday cell
302 mm			Quartz Faraday cell
296 mm			Quartz Faraday cell
280 mm			Quartz Faraday cell
254 mm			Quartz Faraday cell

Specifications


Principle	Automatic digital polarimeter with symmetric angular oscillation using the optical-null balance method
Light source	Tungsten-Halogen lamp (WI), Sodium lamp (Na), Mercury lamp (Hg) (Up to two light sources can be installed.)
Modulator	Faraday cell
Wavelength	880, 633, 589, 578, 546, 436, 405, 365, 334, 325, 313, 302, 296, 280, 254 nm
Aperture	1.8, 3 and 8 mm diameters
Angular range	±90°
Response speed	6°/sec
Measurement accuracy	±0.002° (up to 1°) ±0.2% (larger than 1°)
Repeatability	0.002°
Resolution	0.0001°
Integration time	1 - 100 sec
Detector	Photomultiplier tube (1P28-01) (R928 - option)
Readout modes	Optical rotation, optical specific rotation, concentration, sugar scale (Z), Brix purity, optical purity
Temperature measurement range	0~40°C (minimum display temp.: 0.1°C)
Dimensions	653 (W) x 249 (D) x 364 (H) mm
Weight	Approx. 30 kg
Power requirement	AC100~240V ±10%, 50 or 60 Hz, 80 - 300 W

Optional accessories



PTC-203 Peltier cell holder

The PTC-203 Peltier sample cell holder (air cooled) is designed to accept both rectangular and cylindrical cells with a temperature accuracy of ±0.1°C.

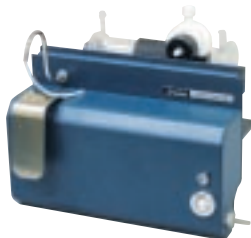

SHP-201W Water jacket sipper

Two types of sample sippers, the SHP-201P Peltier temperature control system and the SHP-201W water-jacketed type for temperature control, are available for quick measurements of multiple samples.



PT-31 Peltier thermostatic bath

The PT-31 is a compact Peltier thermostatic bath whose dimensions are only 140 (L) x 80 (W) x 210 (H) mm. The temperature range is 8 to 40°C with a temperature accuracy of ±0.2°C.



● Specifications are subject to change without notice.

Tecnología Aplicada Internacional

Tel (506) 297 1011/ Fax (506) 297 2231


Apartado Postal 1581-2150

San José, Costa Rica

América Central

www.tecnologia-aplicada.com

info@tecnologia-aplicada.com

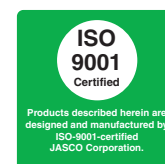

Serving the global marketplace
with Analytical Instrumentation

JASCO Corporation

2967-5, Ishikawa-cho, Hachioji, Tokyo 192-8537, Japan
http://www.jasco.co.jp